

OLYMPUS[®]

Your Vision, Our Future

Endo**Therapy**

Single Use Electrosurgical Snares SD-400

SnareMaster Plus

Cold Snare Polypectomy



SnareMaster Plus

Specially designed for **cold** snare polypectomy

Newly designed hexagonal shape

- The hexagonal shape allows for precise manipulation and control when capturing polyps, whilst maintaining an appropriate margin of normal tissue required for cold snare polypectomy



Durability

- Ability to retain loop configuration (shape) after multiple cuts

For cold and hot polypectomy

- Equipped with additional connector for hot polypectomy if required



Thin wire

- 0.3 mm thin wire with adequate stiffness facilitates clean, sharp cold cutting




Smooth operation

- Special coating of the wire facilitates smooth operation
- Less force is required to cut the lesion sharply and quickly

Size line up (10mm, 15mm)

- Appropriate size line up to resect small polyps by both cold and hot polypectomy

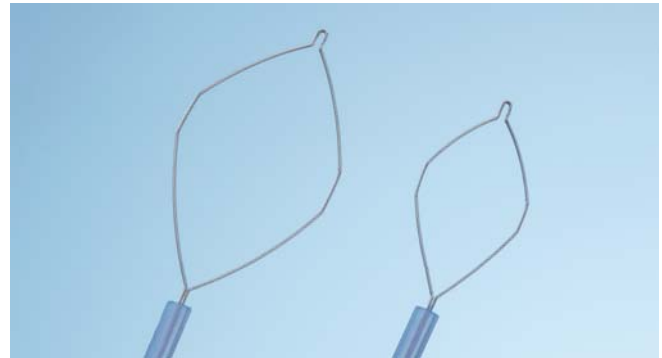

Specifications

Model	Loop design	Working length	Loop diameter	Maximum insertion portion diameter	Diameter of snare loop wire
SD-400U-10	Hexagon	2300 mm	10 mm	2.6 mm	0.3 mm
SD-400U-15			15 mm		

SD-400U-10/15



Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.

OLYMPUS[®]

OLYMPUS AUSTRALIA PTY LTD
3 Acacia Place, Notting Hill VIC 3168, Australia
Customer Service: 1300 132 992 | www.olympusaustralia.com.au

OLYMPUS NEW ZEALAND LIMITED
28 Corinthian Drive, Albany, Auckland, New Zealand
Customer Service: 0508 659 6787 | www.olympus.co.nz